




HYPERTHERMIA PUMP™

Safely Deliver Hyperthermic Therapy

SAVING LIVES. TOGETHER.



INTRAPERITONEAL
& INTRATHORACIC
HYPERTHERMIA

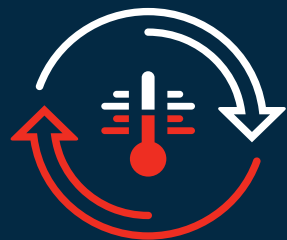

HIGHLY EFFICIENT



USER-FRIENDLY



COMPACT DESIGN

HYPERTHERMIA PUMP™

Hyperthermia Pump™ is used in therapeutic procedures to raise the temperature of the thoracic or peritoneal cavity to the desired target temperature by continuously lavaging the cavity with warmed sterile solution. The warm solution is pumped into a body cavity, withdrawn, reheated, and recirculated back to the body. To ensure safe operation, the system monitors line pressure and air in the fluid path while two sterile temperature probes located within the sterile field and two embedded temperature probes monitor patient cavity, outflow solution, and return solution temperature.

PRECISE TEMPERATURE CONTROL

Adjustable fluid temperature control range between 37°C and 48°C in 0.1°C increments.

ELECTROMAGNETIC INDUCTION HEATING

Dry, single pass warming provides consistent fluid temperature and enables near instantaneous warming to target temperatures with zero warmup time.

COMPACT AND LIGHTWEIGHT

Designed with space in mind, the system conveniently mounts to any standard IV pole and is up to 84% lighter (28.5 vs. 175 lb) and 90% more compact (0.7 vs. 6.8 ft³) than comparable systems.

HIGH-SPEED PERISTALTIC PUMP

Precise control of fluid delivery with flow rates ranging from 2.5 to 1000 mL/min.

VACUUM ASSISTED OUTFLOW

Incorporates vacuum regulator to maximize flow equilibrium and enable fine tuning of fluid return rate.

AUTOMATIC PRESSURE REGULATION


Line pressure monitoring stops pump if line becomes obstructed to prevent leaks.

QUICK SETUP

Disposables incorporate pre-assembled cannulas with quick release connections and embedded temperature probes.

EASY TO USE

High contrast, clutter-free display provides quick access to device operations and system status.

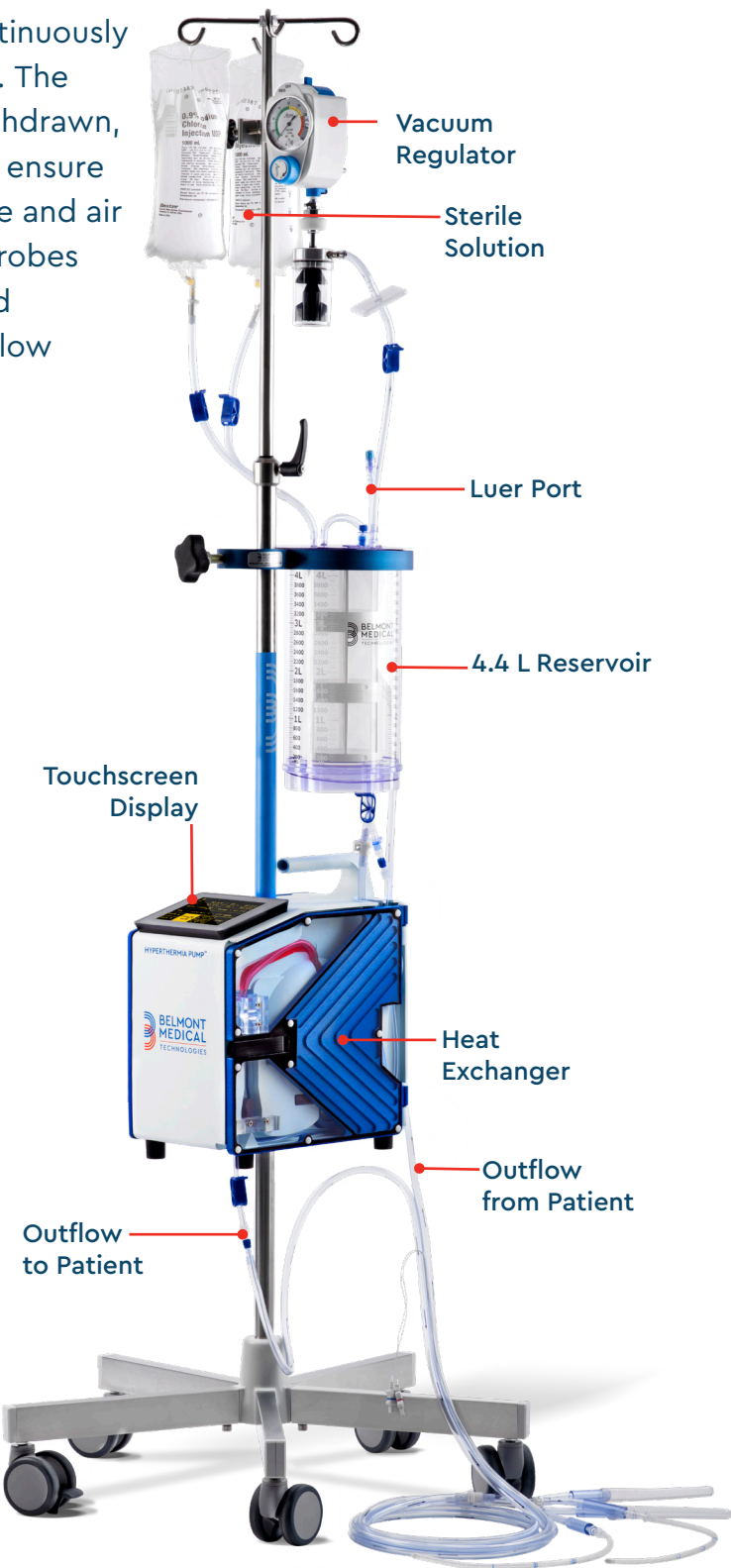

INTUITIVE TOUCH-SCREEN DISPLAY


T1 = 42.0°C	T2 = 42.2°C
T3 = 42.5°C	T4 = 42.6°C
RATE = 1000 $\frac{\text{mL}}{\text{min}}$	T _{pump} = 42.7°C
VOL = 2173 mL	P = 172 mmHg
RATE ▲ 1000 $\frac{\text{mL}}{\text{min}}$	TARGET ▲ 42.5°C
RATE ▼	TARGET ▼ 42.5°C
	STOP

REAL-TIME DISPLAY OF:

- Flow Rate
- Line Pressure
- Target Fluid Temperature
- Total Volume
- Fluid Inflow Temperature
- Fluid Outflow Temperature
- Two Additional Temperature Sensors

100% ALUMINUM-FREE DISPOSABLES

Medical grade stainless steel heat exchanger eliminates risk of aluminum exposure.

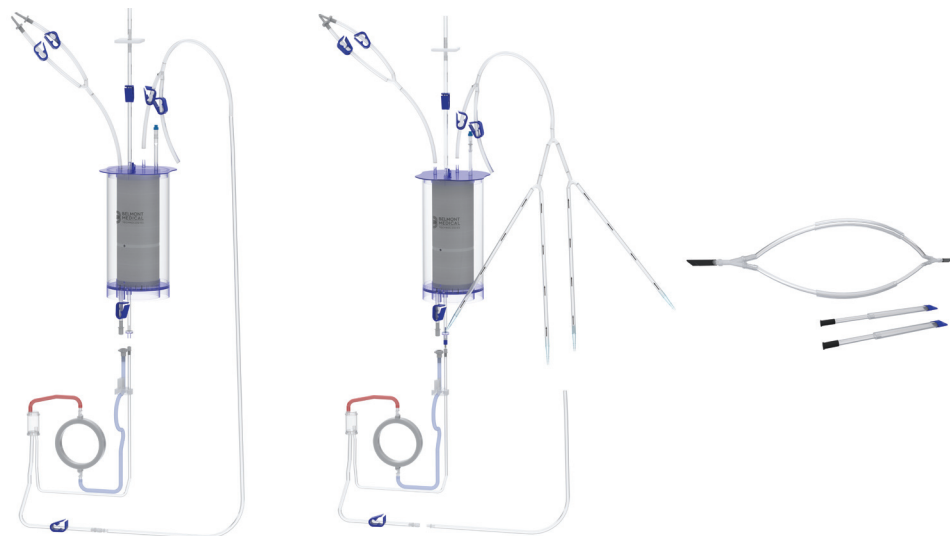


PROCEDURE PACK

The Procedure Pack is designed to save time and money by simplifying complicated tubing setups into a single, comprehensive and easy to install kit.

NO CUTTING, NO CLAMPS – JUST CLICK AND CONNECT.

- Quick-connect luer connections allow for easy setup. There is no need to cut and reconnect tubing.
- Inflow tubing temperature probe reads the inflow temperature before bifurcations.
- Outflow tubing temperature probe automatically averages the temperature of all returned fluids.
- 72" temperature probes allow for connection outside of the sterile fluid field.

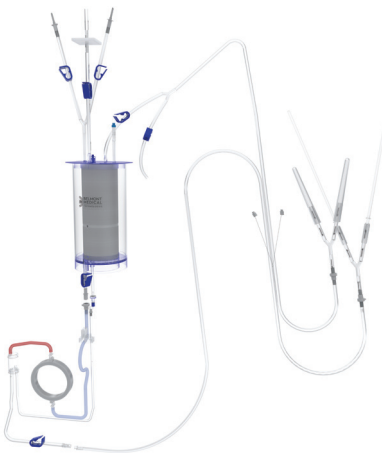

ORDERING INFORMATION

HYPERTHERMIA PUMP

- US: Part #902-00001 (120V)
- International: Part #902-00001A (230V)

CANNULA CONNECTOR

Part #902-00036 (12/box)

STRAIGHT LINE SET

Part #902-00037 (1 kit)

BIFURCATED LINE SET

Part #902-00038 (1 kit)

PROCEDURE KIT

Part #902-00045

(1) 4.4 Liter Reservoir (902-00034P)

(1) Heat Exchanger (902-00006P)

(1) Procedure Pack (902-00048P) contains:

Double lumen patient line with (2) temperature probes*

- (2) Extra temperature probes
- (2) Sumps
- (2) Via-Guard®

*Each temperature probe is 72 inches (1.83 meters) in length, allowing connection outside the sterile field.

ALL DISPOSABLES ARE 100% ALUMINUM-FREE

HARDWARE SPECIFICATIONS

DIMENSIONS (INCLUDING FEET AND HANDLES)

12.4 in. x 7.5 in. x 14.8 in. (31.5 cm. x 19.1 cm. x 37.6 cm.)

WEIGHT

28.75 lb (13.04 kg)

PORTABILITY

Quick-on IV pole mount/handle at top of unit

POWER

115–120 volt or 220/240 volt AC input

BATTERY

Internal, rechargeable for transport

Operates all functions without heating

FLUID PUMP

Roller type peristaltic pump

Flow rate 2.5-1000 mL/min

HEATER

Dry, electromagnetic induction (water-free)

Adjustable between 37–48°C in 0.1°C increments

Maximum applied power: 1440 watts

SENSORS

Up to 4 temperature sensors (2 patient cavity, 1 outflow,

1 return), 2 ultrasonic air detectors, pressure transducer,

2 infrared temperature sensors, pump infusion rate

sensor, open door detector, and valve activation sensors

CONTROL COMPUTER

Built-in, controls all functions

Redundant hardware override protection in case

of computer failure

DISPLAY

LCD type with waterproof touch pad

• External Dimensions

5.85 in. x 4.12 in.

148.1 mm x 104.7 mm

• Active Area

4.53 in. x 3.40 in.

115.1 mm x 86.3 mm

For more information about the products we offer,
contact us at sales@belmontmedtech.com.



PATIENT
WARMING



TEMPERATURE
MANAGEMENT



FLUID
MANAGEMENT



780 BOSTON ROAD
BILLERICA, MA 01821, USA
USA: 866.663.0212
WORLDWIDE: +1 978.663.0212
WWW.BELMONTMEDTECH.COM



Hyperthermia Pump™ has been evaluated with respect to electrical shock, fire and mechanical hazards in accordance with applicable standards.

Belmont Instrument LLC
Registered to ISO 13485
File Number: A18173

Belmont Medical Technologies™ is a trademark of Belmont Instrument, LLC
2023 © Belmont Instrument, LLC | 701-00014 Rev. 001