



# Quick-Core<sup>®</sup> Biopsy Needle and Needle Set

Minimally Invasive Soft Tissue Biopsy

## What is Quick-Core<sup>®</sup> Biopsy Needle?

The Quick-Core<sup>®</sup> Biopsy Needle is a device used in minimally invasive soft-tissue biopsy procedures that can be performed in a specialists office or in a hospital setting.

The Quick-Core<sup>®</sup> Biopsy Needle is indicated for soft-tissue biopsy.

## How Quick-Core<sup>®</sup> Biopsy Needle and Needle Set Works

The soft tissue biopsy procedure can be done under local anesthetic. The procedure usually takes less than an hour to complete.

The Quick-Core<sup>®</sup> set includes a stainless steel trocar.

## Why use Quick-Core<sup>®</sup> Biopsy Needle versus Competitors?

- Smart design reduces moving components allowing multiple uses without breakage
- Sharp bevel tip with two additional facets facilitates easy placement into tissue
- Rigid spring reduces vibration which allows better cutting of specimen preserving tissue architecture
- Small, ergonomic and lightweight handle limits needle movement
- No device recalls and a legacy of quality



Quick-Core® / Needle

# Quick-Core® Biopsy Needle

## Soft-tissue Biopsy Needle

Used for soft-tissue biopsy

- The lightweight and compact Quick-Core® Biopsy Needle **facilitates one-handed operation** during procedures
- The **precise, rapid-firing mechanism** provides an intact core tissue sample
- The beveled stylet **permits easy penetration** into the lesion with less trauma to surrounding tissues and organs
- The etched tip **enhances visualization** under ultrasound


Supplied sterile in peel-open packages. Intended for one-time use only.



ORDER INFORMATION

**Quick-Core® Biopsy Needle**

SET ORDER NUMBER	OUTER CANNULA GAUGE • LENGTH	GLOBAL PROD. NO.
<b>10 MM THROW LENGTH</b>		
QC-16-6.0-10T	16G • 6 cm	G28280
QC-16-9.0-10T	16G • 9 cm	G28282
QC-18-6.0-10T	18G • 6 cm	G28283
QC-18-9.0-10T	18G • 9 cm	G27009
QC-18-15.0-10T	18G • 15 cm	G07819
QC-20-6.0-10T	20G • 6 cm	G27012
QC-20-15.0-10T	20G • 15 cm	G07822
<b>20 MM THROW LENGTH</b>		
QC-14-9.0-20T	14G • 9 cm	G27006
QC-14-15.0-20T	14G • 15 cm	G07824
QC-16-6.0-20T	16G • 6 cm	G28281
QC-16-9.0-20T	16G • 9 cm	G27007
QC-16-15.0-20T	16G • 15 cm	G07794
QC-16-20.0-20T	16G • 20 cm	G07827
QC-18-6.0-20T	18G • 6 cm	G27008
QC-18-9.0-20T	18G • 9 cm	G27010
QC-18-15.0-20T	18G • 15 cm	G07778
QC-18-20.0-20T	18G • 20 cm	G07821
QC-20-9.0-20T	20G • 9 cm	G27015
QC-20-15.0-20T	20G • 15 cm	G07780

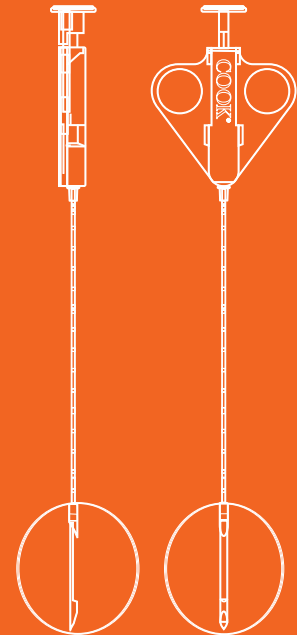

THROW LENGTH

Not all parts numbers shown may be approved for sale in all regulatory jurisdictions. Parts and configurations are subject to change without notification. Please consult with your Account Representative for details.

Quick-Core® / Set

## Quick-Core® Biopsy Needle Set

### Soft-tissue Biopsy Needle Set

Used for soft-tissue biopsy

- Coaxial outer needle permits multiple biopsies with the Quick-Core® needle through a **single access site**
- The lightweight and compact Quick-Core® Biopsy Needle facilitates **one-handed operation** during procedures
- The **precise, rapid-firing mechanism** provides an intact core tissue sample
- The EchoTip® surface of the outer coaxial needle and the etched tip of the cutting cannula **enhance visibility** under ultrasound


Supplied sterile in peel-open packages. Intended for one-time use only.



ORDER INFORMATION

**Quick-Core® Biopsy Needle Set**

SET ORDER NUMBER	BIOPSY NEEDLE GAUGE • LENGTH	COAXIAL NEEDLE GAUGE • LENGTH	GLOBAL PROD. NO.
<b>10 MM THROW LENGTH</b>			
QCS-18-9.0-10T	18G • 9 cm	16G • 6 cm	G27016
QCS-18-15.0-10T	18G • 15 cm	16G • 12 cm	G08728
<b>ULTRA-THIN WALL</b>			
QCS-20-9.0-10T	20G • 9 cm	19G • 6.5 cm	G27018
QCS-20-15.0-10T	20G • 15 cm	19G • 12.5 cm	G08727
<b>20 MM THROW LENGTH</b>			
QCS-14-9.0-20T	14G • 9 cm	11G • 5.4 cm	G28290
QCS-18-9.0-20T	18G • 9 cm	16G • 5 cm	G27017
QCS-18-15.0-20T	18G • 15 cm	16G • 11 cm	G08790
QCS-18-20.0-20T	18G • 20 cm	16G • 16 cm	G08789
<b>EXTRA-THIN WALL</b>			
QCS-16-9.0-20T	16G • 9 cm	14G • 5.6 cm	G28292
QCS-16-15.0-20T	16G • 15 cm	14G • 11.6 cm	G10149
QCS-20-9.0-20T	20G • 9 cm	19G • 5.5 cm	G27019
QCS-20-15.0-20T	20G • 15 cm	19G • 11.5 cm	G08788
QCS-20-20.0-20T	20G • 20 cm	19G • 16.5 cm	G08787

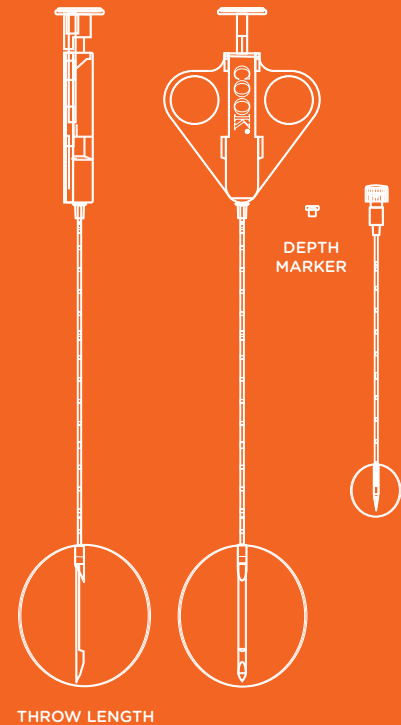

Not all parts numbers shown may be approved for sale in all regulatory jurisdictions. Parts and configurations are subject to change without notification. Please consult with your Account Representative for details.